



Technical resources

Laboratory sand grinder

Model: SM

<https://www.tencanballmills.com/en/products/grinding-series/other-grinding-equipment/laboratory-sand-grinder.html>

09.10.2026 / TENCAN

Laboratory sand grinder


Model: SM

Multi-functional lab sand grinder for crushing, dispersing, and emulsifying metals, non-metals, organics, and herbal powders. Ideal for R&D and process simulation.

Technical specifications

Agitator	2 grinding discs
Agitator speed	1400~3500r/min
Ground particle size	Can reach 1 μm

These specifications describe the series. Confirm the selected model, power supply and configuration before ordering.

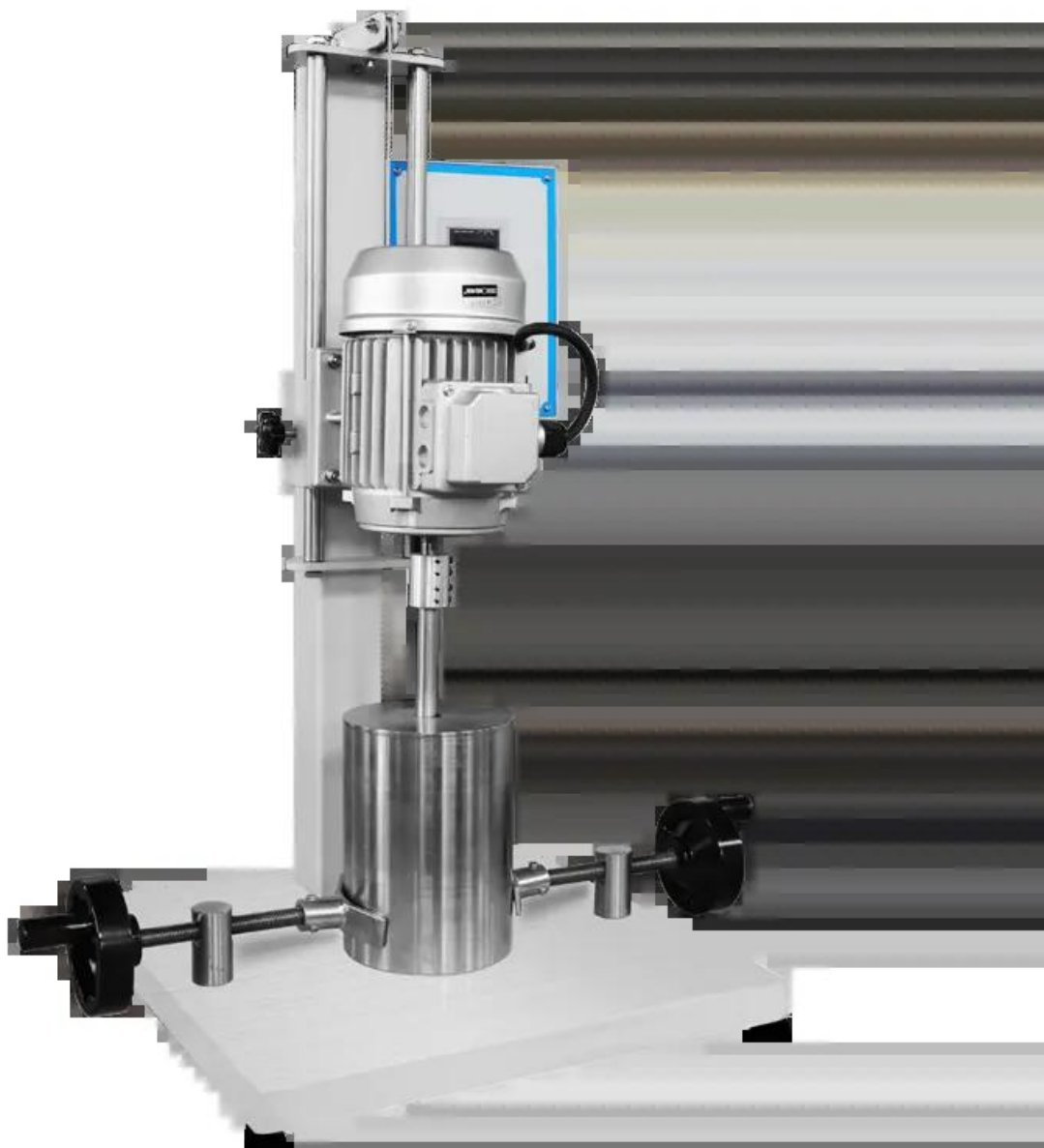

The laboratory sand mill is also called a dispersing sand grinding machine or a dispersing sand grinding dual-purpose machine. It is designed for crushing, grinding, dispersing, emulsifying metal, non-metallic, organic, Chinese herbal medicine and other powders. It is especially suitable for laboratory research. It can well achieve various process parameter requirements and simulate various indicators in production. At the same time, due to its advantages of small batch, low power consumption and low price, it is an optional equipment for schools, research units and companies to conduct research on crushing processes, new materials and coatings.

Experimental sand grinding and dispersing machines are mainly used in paint, ink, coating chemical industry and electronic ceramics, structural ceramics, magnetic materials, lithium cobalt oxide, lithium manganate, catalysts, phosphors, long afterglow luminescent powder, rare earth polishing powder, electronic glass powder, fuel cells, ceramic capacitors, oxygen Zinc varistors, piezoelectric ceramics, nanomaterials, disc ceramic capacitors, MLCC, thermistors (PTC, NTC), ZnO varistors,

arrester valves, strontium titanate ring varistors, ceramic filters, dielectric ceramics, piezoelectric transducers, piezoelectric transformers, chip resistors, thick film circuits, pyropotentiometers, alumina ceramics, zirconia ceramics, Production and research industry of phosphors, zinc oxide powder, cobalt oxide powder, Ni-Zn ferrite, Mn-Zn ferrite and other products.

The working principle of the experimental sand grinding and dispersing machine is to use the abrasive and the sample to tumble at high speed in the grinding tank to produce strong shearing, impact and rolling of the material to achieve the purpose of crushing, grinding, dispersing and emulsifying the material.







Technical parameters

Model specifications	unit
Grinding jar volume	0.5□1□2□3□5□10L
Grinding jar lining	Stainless steel, alumina, zirconia, etc.
mixer	2 sanding discs
Stirring speed	1400, 2800, 3500r/min or adjustable speed
Mosuke	Suitable for grinding media such as corundum and zirconia below Φ10
Motor power	120-1500W
Crushing particle size	Up to 1μm

1. Simple structure and easy operation;
2. The mixer adopts a quick-change chuck, which is easy to install, remove and replace. ;
3. Use a sanding barrel with a cooling jacket, which can be cooled while grinding. ;

4. Replace the agitator with a dispersion plate to disperse the logistics. ;
5. The material can be ground at different speeds according to process requirements.

Contact

Changsha Tianchuang Powder Technology Co., Ltd. / Changsha, Hunan, China

Email: lichang0516@gmail.com

Phone / WhatsApp: +86 731 84027560 / WhatsApp +86 19021571516

<https://www.tencanballmills.com/en/products/grinding-series/other-grinding-equipment/laboratory-s-and-grinder.html>

Before operation, review the selected model's instructions and the material requirements.